


At Last... A Wireless Nurse Call System Approved for Skilled Care!

VisionLink™ Wireless Nurse Call System Flexible, Reliable, Affordable

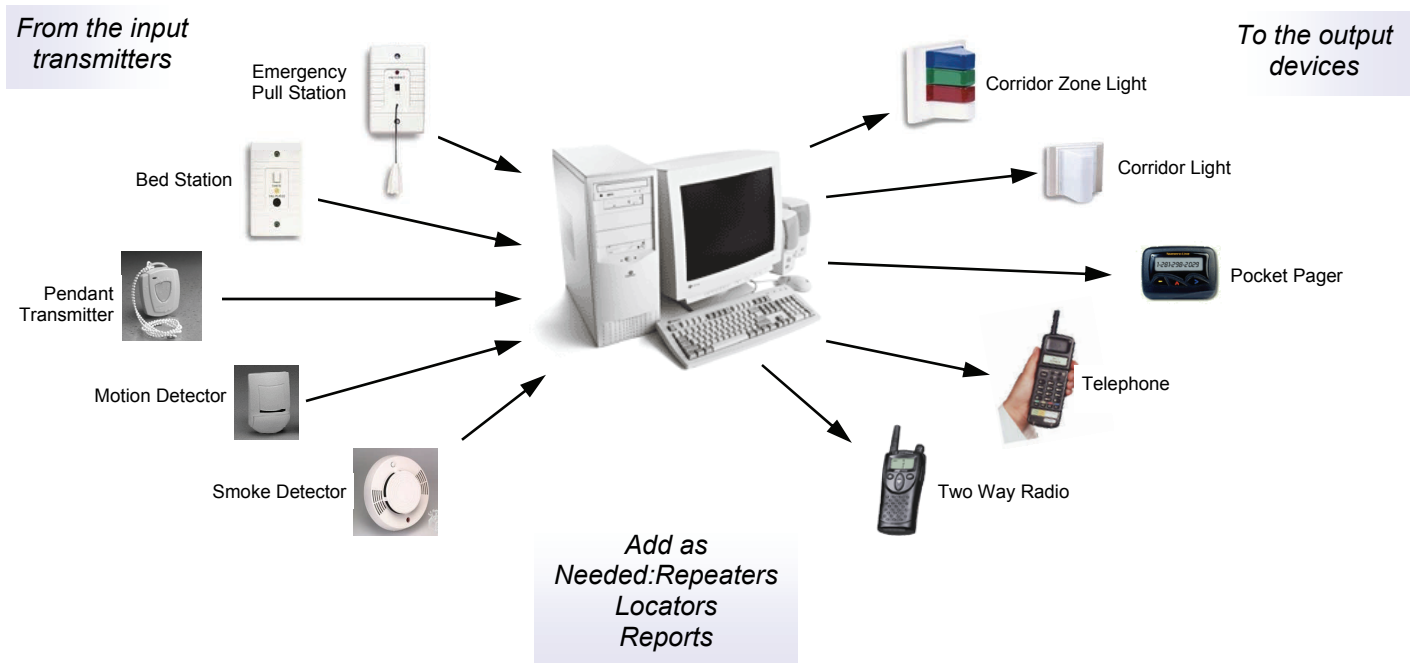

Start Small, Think Big!

Vision Link utilizes the Windows 2000® platform combined with 900 MHz spread spectrum wireless technology to link your residents to your staff. Vision Link gives you the ability to start small and think big. Go wireless when adding a new wing or new building to save on conduit and wiring costs. For skilled care, add a single two wire loop through the corridors to provide dome lights. For a multi-building campus, add repeaters to send assistance calls to the nearest console by day and a centralized console at night.

Flexible Software

Patient information is a standard Vision Link feature along with a call history log that records every resident request for assistance and staff response times. Also standard are Call Analysis Reports that provide the tools to organize and analyze the data. Combination Emergency Pull/Check-In stations allow residents to summon help if needed and notify staff they are up and active each day. A scheduler reminds staff to give medications, lock doors or perform other tasks. Bed turn software notifies CNA's to turn their patients according to schedule.

Staff Mobility

Staff members are no longer bound to the console. Caregivers can use pocket pagers, two way radios or wireless phones to receive alarm calls. Vision Link features text messaging. Not just patient alarms but any message typed on a console keyboard can be sent to any pager or pager group eliminating institutional overhead paging. Staff can be text messaged that calls are holding or other information. Vision Link can also provide notification from wandering patient systems, exit doors, smoke detectors, motion detectors, bed exit pads, equipment status, whatever your facility requires.

The 900MHz Advantage

Because the FCC allows higher transmit power for 900 MHz spread spectrum systems, Vision Link provides superior in-building range and up to 10 times the range of 300 MHz products. It is also fully supervised and transmits every 10 seconds to indicate low battery or inactive status. In contrast, it can take some wireless systems up to 24 hours to report a malfunctioning transmitter.

Improve safety for residents and staff with a *VisionLink™* Wireless Nurse Call System

EQUIPMENT OPTIONS

- Wireless Patient, Code and Duty Stations
- Wireless Emergency Pull Stations
- Wireless Pendant and Wrist Transmitters
- Corridor and Zone Lamps (white or two color)
- Alpha-numeric pocket page alarm notification
- Central station alarm monitoring interface
- Telephone and text message interface
- Networked consoles can select response areas
- Door/Window/Bed exit notification
- Wandering patient notification
- Smoke/motion/glass breakage detection
- Resident check-in status
- Repeaters for campus coverage
- PC or standard console
- locators to pinpoint alarm source

Date/Time	State	Type	Room	Bed	Occupant
2/1/2002 4:15:04 PM	Call	Smoke Det.	302		
2/1/2002 4:15:08 PM	Call	Smoke Det.	Lounge		
2/1/2002 4:16:03 PM	Call	Smoke Det.	103	1	Angela Bennett
2/1/2002 4:15:04 PM	Call	Code Blue	205		MARY-JONES
2/1/2002 4:15:08 PM	Call	Code Blue	202	2	Marisa Shelton
2/1/2002 4:16:03 PM	Call	Code Blue	101	1	Chuck Smith
2/1/2002 4:15:04 PM	Call	Pull Station	201	2	Sue Davis
2/1/2002 4:15:04 PM	Call	Pendant	202		
2/1/2002 4:15:08 PM	Call	Pull Station	201		
2/1/2002 4:16:03 PM	Call	Pull Station	203		
2/1/2002 4:15:04 PM	Call	Pendant	102	1	Gordon Hutson
2/1/2002 4:15:08 PM	Call	Pendant	105	1	Holly Smoke
2/1/2002 4:15:04 PM	Call	Bed Station	303	1	Wade Shepard
2/1/2002 4:15:08 PM	Call	Bed Station	104	1	Sally Petruska
2/1/2002 4:16:03 PM	Call	Bed Station	101	1	Chuck Smith

Call Today

Our experienced, knowledgeable technical staff would be happy to answer your questions or give you ideas on how to put wireless to work for your facility.



STANDARD FEATURES

- Windows 2000® operating system
- 900 MHz spread spectrum technology
- Fully supervised (self-testing)
- Alarm Calls reset at the source
- Low battery alert
- Password protection
- Visual/audible alarm prioritization
- User replaceable long-life battery

MORE FEATURES

- Pocket pager alarm notification
- Two way radio alarm announcement
- Programmable page escalation
- Specialized management reports
- Passive/active check-in status
- Quality assurance verification
- Timed egress alarm notification
- Text messaging
- Staff presence indication
- Bed turn scheduler

Bed	Occupant	Location
101	Chuck Smith	North Wing



Systemstechnologies: Transforming Technology Into Solutions

For more information, call or write:

NurseCall Northwest
Phone: 206-686-7111
Fax: 206-686-7111

Address: 6523 California Street SW, PMB 107
Seattle, WA 98136
Email: sales@nursecallnw.com



Systems Technologies reserves the right to change specifications without notice. Some features may require additional hardware or software. © 2005 Systems Technologies, Inc. All rights reserved. Patents pending worldwide. All products or services mentioned herein are the trademarks of Systems Technologies, Inc., or their respective owners.



UL1069 Listed